

novar *FREEFORM LENS DESIGNER*


INDEX

TECHNOLOGIES

- 03 CONTINUUM DESIGN TECHNOLOGY
- 04 4K DEFINITION SURFACE
- 05 ADAPTATIVE FOCUS TECHNOLOGY
- 06 WEAR FIT
- 07 SMART BLENDING
- 08 AUTOMATIC CORRIDOR
- 09 SMART FIT
- 10 RAYTRACING

DESIGNS

- Monofocals*
 - 11 SINGLE VISION
 - 12 GEO
 - 13 GEO SMART
- Ocupacionals
& Bifocals*
 - 14 OFFICE
 - 15 KRIPTOK BLENDED
 - 16 ULTEX BLENDED
 - 17 FREELINE
- Multipurpose
Progressives*
 - 19 FIRST
 - 20 ECOLINE
 - 21 PRECISA
 - 22 PRECISA SHORT
 - 23 EVOLUTION
 - 24 EVOLUTION SHORT
 - 25 eLIFE
- Specialty
Lens*
 - 26 RELAX
 - 27 DRIVE
 - 28 OUTDOOR
 - 29 INDOOR
 - 30 MONOVISION
 - 31 SENSITIVE
 - 32 KIDCARE
- 33 SPECIFICATIONS & CHARACTERISTICS OF THE DESIGNS
- 34 THICKNESS TECHNOLOGY
- 35 BLENDING TECHNOLOGY
- 37 SOFTWARE DESIGNER
- 42 MARKETING
- 44 DIGITAL CENTER

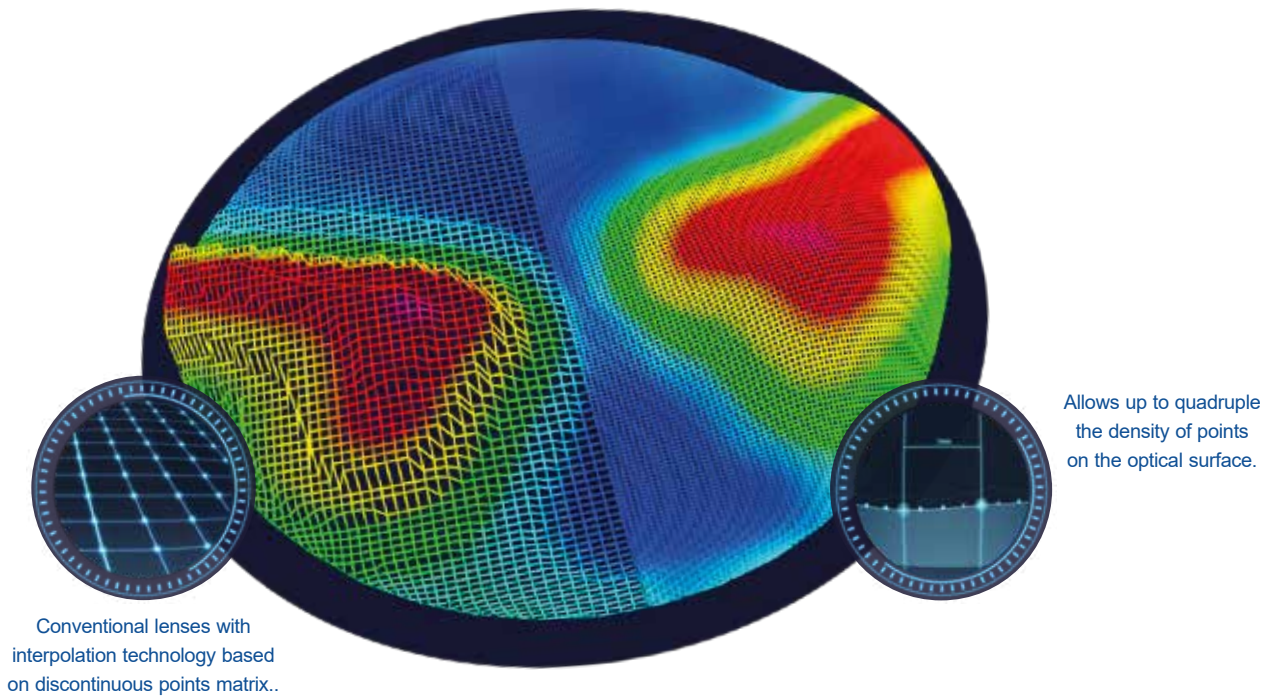

CONTINUUM DESIGN TECHNOLOGY

CDT is an innovative design technology based on a modern mathematical model that allows improving the manufacture of ophthalmic lenses.

With this breakthrough technology, NOVAR lenses have a starring role in the new technological process within the global market.

While other technologies are only based on defining optical conditions within a certain set of points on the lens surface and then, interpolate (b-splines), CONTINUUM DESIGN TECHNOLOGY allows to take control over the whole lens surface to send greater matrices of sagittas to the generator. Consequently, a higher optical resolution to optimize visual fields is achieved.



*Sharper & brighter
Images*



4K DEFINITION SURFACE

CONTINUUM DESIGN TECHNOLOGY allows to produce lenses with up to 4 times more resolution than other Freeform technology lenses.

As pixels in U-HD digital screens, a greater resolution on the optical surface turns into real benefits linked to greater visual field amplitude and image resolution.




ADAPTATIVE FOCUS TECHNOLOGY

used in our new elife progressive lens

The human being has been part of radical changes throughout its evolutionary history. Thus, we went from reading on carved stones to reading on modern digital devices at an overwhelming speed. Speed that involved a substantial change in the way we see things.

As digital media such as smartphones, tablets, etc grew exponentially, we began to adopt focus habits that challenge the adaptability of our eyes. This increased the risk of suffering from visual fatigue, neck pain, headache or even led us to Computer vision Syndrome (CVS). Consequently, we developed Adaptive Focus technology, in which the intermediate and near visual field, most commonly used for reading digital and printed media, has been optimized.

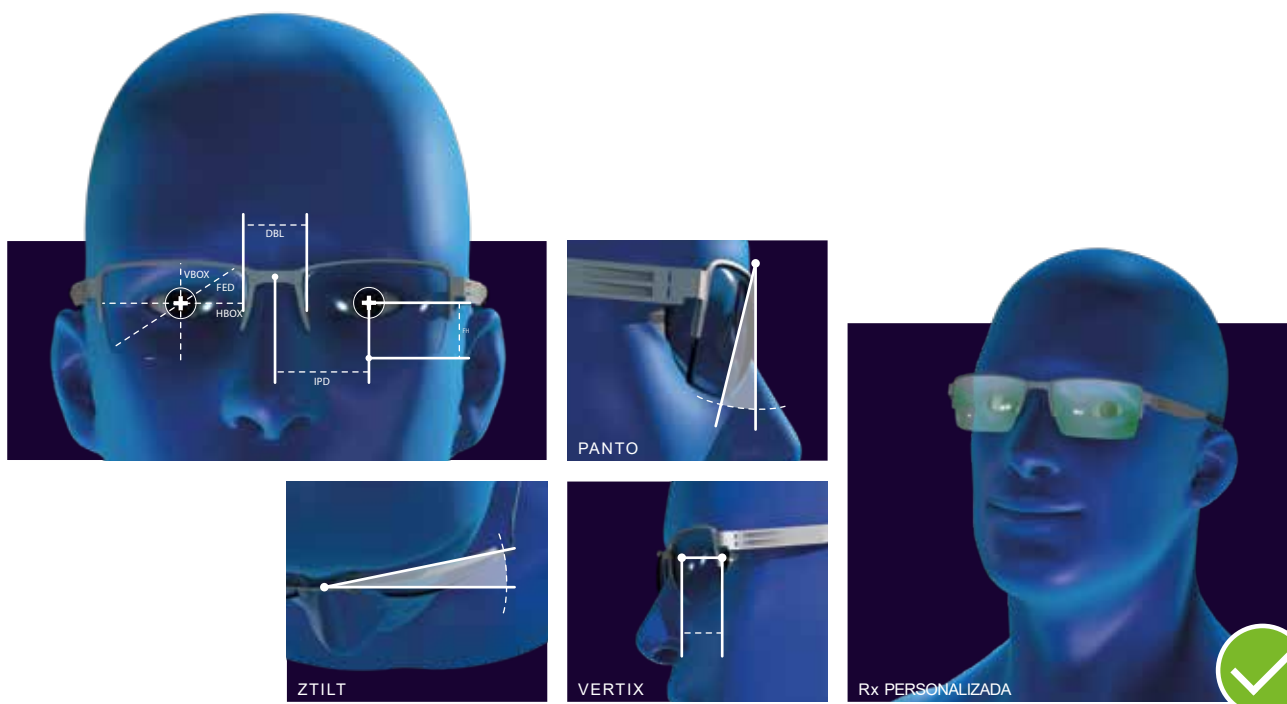

WEAR FIT

The new era led us to be more and more demanding and get the best out of the things surrounding us. However, many times we are forced to adapt pre-established patterns that do not fit our peculiarities at all. For instance, while choosing a frame, we are acquiring an important but without-customization accessory. Thus, when using the prescribed lenses, the wearer may experiences distortions derived from the position of wear (POW) that were not taking into account and could significantly affect lens performance. This key factor, led us to develop a complex system based on mathematical algorithms which are used in the design and the production process of lenses. It consists in a precise data entry process where all the wearer measurements such as interpupillary distance (IPD), panoramic angle (ZTILT), vertex distance (BVD) and pantoscopic tilt (PANTO) are stored. Consequently, this technology enable us to produce lenses that fit you as a quality custom made suit.





SMART BLENDING

Aesthetics has a starring role in our lives, especially when it is related to health. Producing comfortable lenses that in turn are aesthetic is an arduous task that can only be achieved by means of the use of state-of-the-art technology. Edge thickness (nearsightedness) and central thickness (farsightedness) reduction is achieved by using this technology.

The circular blending process establishes an optical zone around the fitting cross to reduce gradually and uniformly the peripheral thickness.

The Smart blending process, unlike other processes, takes into account lens zones which are mostly used by the visual system. Thus, a variable and Smart blending is applied to reduce every biologically untapped vision zone. As if it were a handcrafted piece, this technology intelligently optimizes thicknesses taking into account the type of frame, the lens prescription and the lens diameter.



AUTOMATIC CORRIDOR

Decision making is a complicated process that requires time to evaluate different options. The lack of time or professional expertise can make this process even harder. Thus, when choosing a frame, for instance, we must take into account its intrinsic characteristics and its interaction with the lenses prescribed for presbyopia. Accordingly, we developed a mathematical algorithm capable of establishing the ideal corridor considering fitting heights, pantoscopic tilt and vertex distance to avoid any problem caused by a wrong decision while choosing the corridor. Thus, maximum comfort and versatility in every vision zone is achieved.



SMART FIT


Accuracy in products manufacturing was and will be a key factor to achieve excellence. This was not only taken into account by old swiss watch brands but also by optical lens manufacturers. In the calculation process, for instance, accuracy led us to a greater thickness reduction resulting in: + OPTICAL QUALITY + AESTHETIC.

Nowadays, thickness optimization is intimately linked to the lens diameter, frame measurements and lens wearer data. SMART FIT, conversely, adds complex calculations that are nourished by the trace shape of the frame. Clearly Speaking, it takes into account the frame geometric shape (All points that form the circumference -TRCFMT) to achieve 40% thickness reduction when finishing the carving process.



RAY TRACING

Raytracing is the standard technique used for lens designs in general and for aspheric and progressive lens designs in particular. It is a technique based on tracing a ray of light through a system by calculating the angle of refraction/reflection at each surface. This enable to optimize its shape according to the required optical properties.

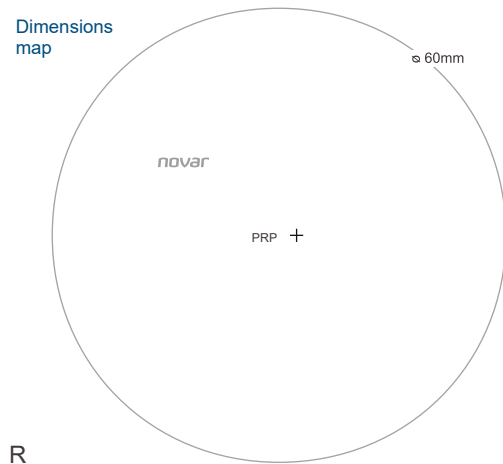Dimensions
map

SINGLE VISION

Spheric/toric monofocal with all the digital quality offered by Freeform.

TECHNICAL DATA



Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

VISION

VISION

Thickness calculation technology:

Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

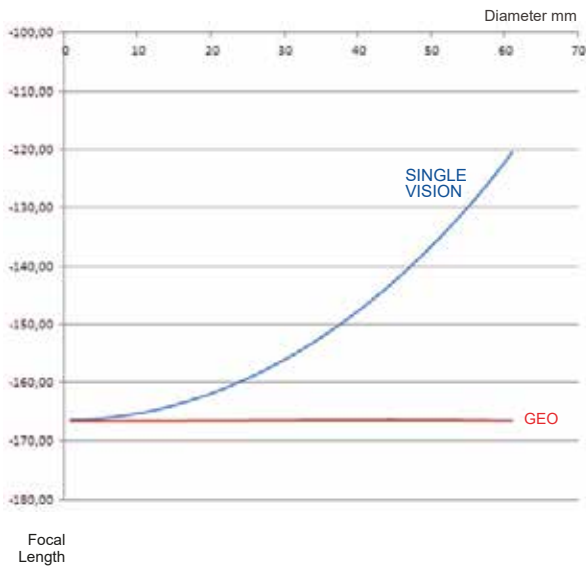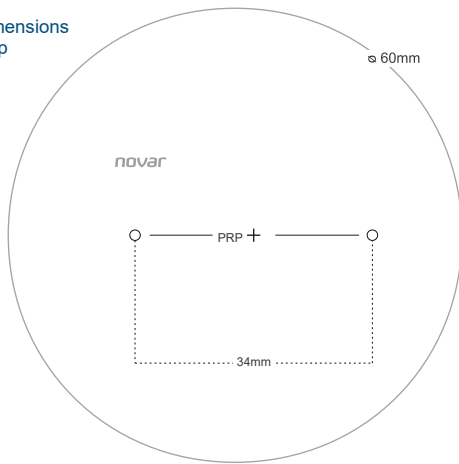Dimensions
map

VISION

VISION

GEO


Aspheric and atoric personalized monofocal which achieves a better visual quality and a higher perception of details by eliminating a great part of spherical aberrations. Consequently, an improvement on visual fields for different gaze directions is obtained.

TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

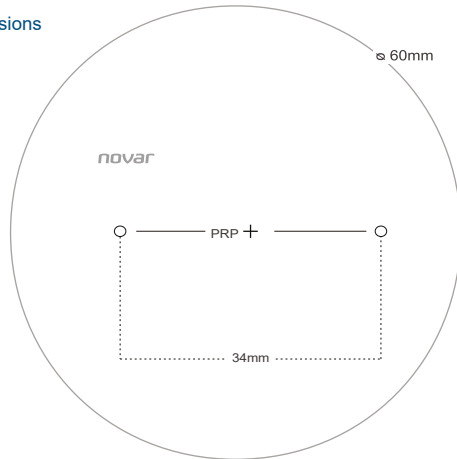Dimensions
map

R

GEO

GEO SMART

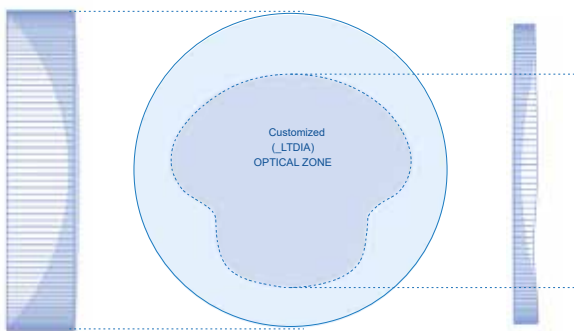Thickness: 8.21 mm
Weight: 22.4 gThickness: 6.23 mm
Weight: 19.8 g

GEO SMART

Aspheric and atoric personalized monofocal which achieves a better visual quality and a higher perception of details by eliminating a great part of spherical aberrations. Untapped visual zones are smartly detected to efficiently improve the thickness of the lens.

TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Maximum diameter	85 mm
Spherical power range	-25 / +25 D
Cylindrical power range	-6 / +6 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positive and Negatives)	Yes

VISION

VISION




Dimensions map

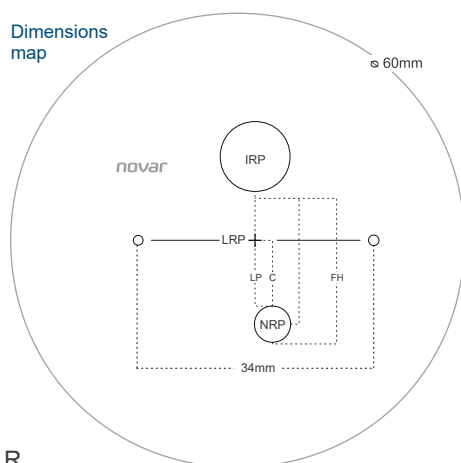


R

Power map

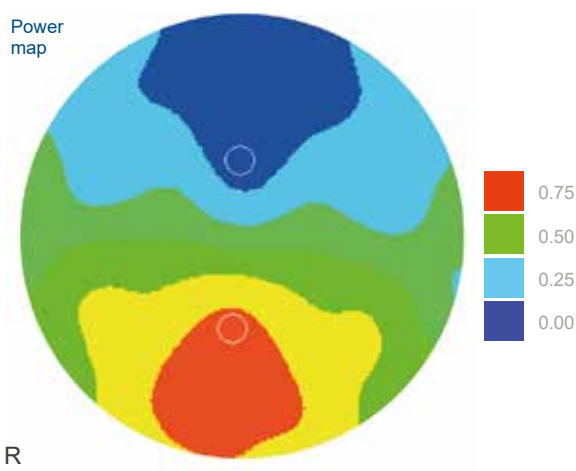

R

Cylinder map

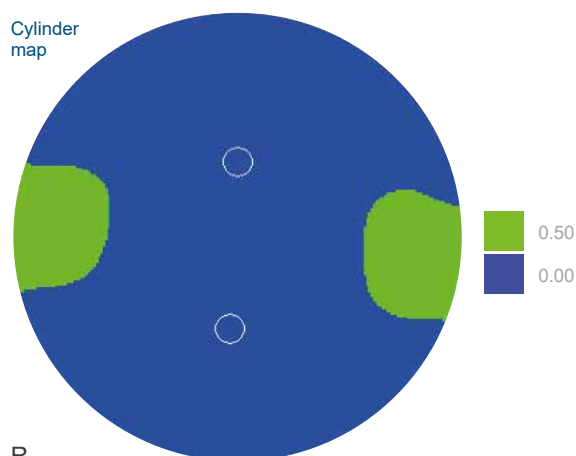

R

VISION



OFFICE

NOVAR Office is a tailor-made lens designed for intermediate and near vision. Suitable for people who need a wider area for near and intermediate vision without peripheral restraints for office work.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Layout reference point (LRP)	Geometric center
Inset	Variable
Minimum VBOX	29 mm
Minimum fitting height (FH)	16 mm
Corridor	24 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.75 / 4.00 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes
Possibility of calculation by degression	0.75, -1.00, -1.25, -1.50, -1.75, -2.00, -2.25
Possibility of calculation by distances	Computer (0,75 mts) Desktop (1,3 mts) Life (1,5 mts) Meeting (2 mts) Room (4 mts)

Thickness calculation technology:

Circular Fit	Yes
Elipitical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

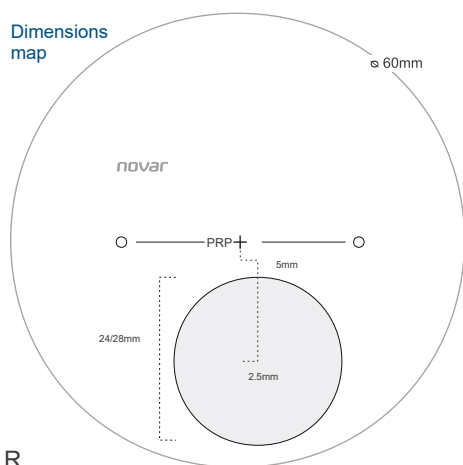

KRIPTOK BLENDED

Specially designed for users who are looking for bifocal designs. Aesthetically improved thanks to its invisible segment + accurate as lens made with Freeform Technology.

TECHNICAL DATA


Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	Geometric center from 0 - 10 mm
Layout reference point (LRP)	Geometric center
Inset	2,5 mm
Segment diameter	24 / 26 / 28 mm
Segment transition	Customized and variable
Minimum fitting height (FH)	14 mm
Vertical displacement	5 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 4.00 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:

Circular Fit	Yes
Elipitical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

VISION



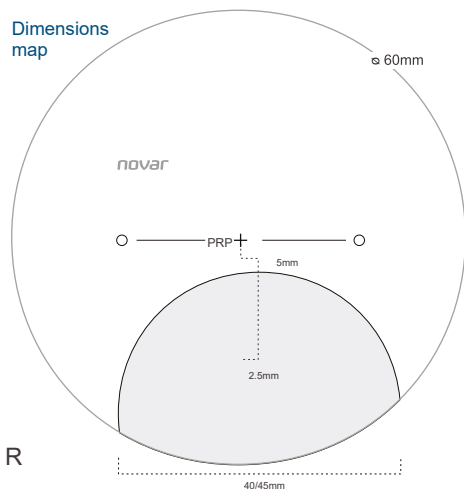

ULTEX BLENDED

Specially designed for users who are looking for bifocal designs. Aesthetically improved thanks to its invisible segment + accurate as lens made with Freeform Technology.

TECHNICAL DATA


Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	Geometric center from 0 - 10 mm
Layout reference point (LRP)	Geometric center
Inset	2,5 mm
Segment diameter	40 / 45 mm
Segment transition	Customized and variable
Minimum fitting height(FH)	14 mm
Vertical displacement	5 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 to 3.00 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:

Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

VISION

FAR	
INTERMEDIATE	
NEAR	

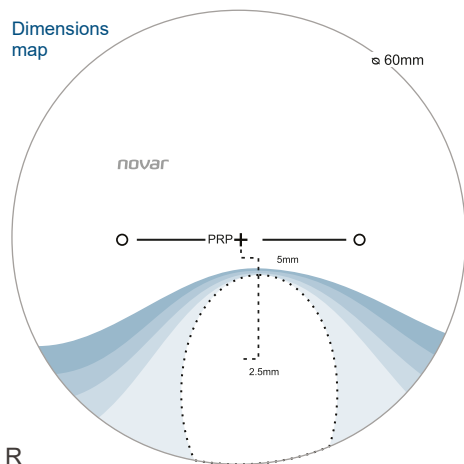

FREELINE

Bifocal with greater aesthetics using the peripheral areas of the lower meridian to improve the invisibility of the segment.

TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	Geometric center from 0 - 10 mm
Layout reference point (LRP)	Geometric center
Inset	2.5 mm
Diameter segment	Aprox. 18mm
Minimum fitting height (FH)	14 mm
Variable decentration	5 mm
Near reference point (NRP)	12 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.00 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:

Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

VISION





Multipurpose Progressive Lenses

STANDARD



FIRST
PROGRESSIVE



ECOLINE
PROGRESSIVE

HIGH



PRECISA
PROGRESSIVE

PREMIUM



EVOLUTION
PROGRESSIVE



eLIFE
PROGRESSIVE

FIRST PROGRESSIVE

For first-time progressive lens wearers. *Distant Clarity Process* is used to maximize far vision breadth and thus achieve a high-performance adaptation. These progressive lenses are ideal for people aged 40-45 years who begin to experience presbyopia.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	Variable
Minimum VBOX	24 mm
Minimum fitting height (FH)	16 - 17 - 18 - 19 - 20 - 21 - 22 mm
Corridor	12 - 13 - 14 - 15 - 16 - 17 - 18 mm
Near reference point (NRP)	14 - 15 - 16 - 17 - 18 - 19 - 20 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.00 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes
Automatic corridor selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

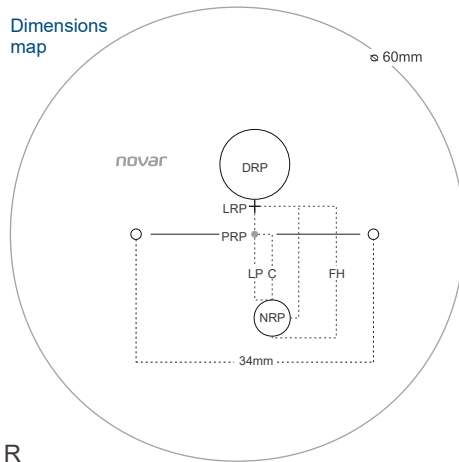


R

Power map

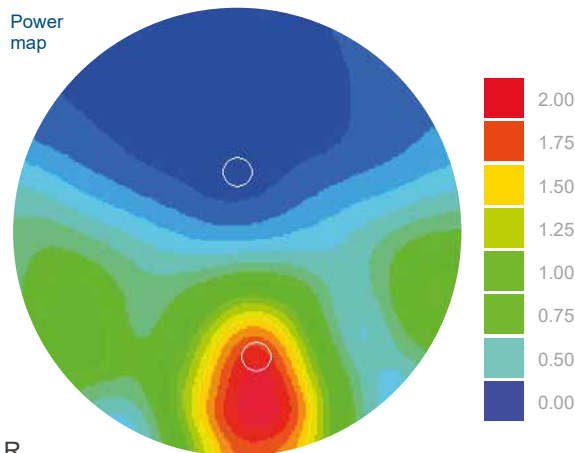

R

Cylinder map

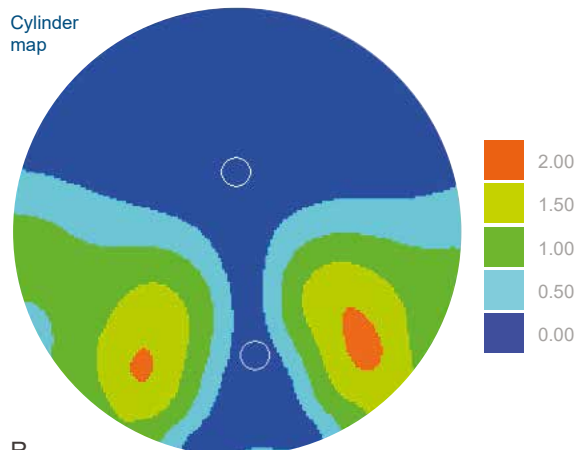

R

VISION



ECOLINE PROGRESSIVE

Progressive lens with excellent price-quality ratio, that provides its users with a great visual performance without compromising quality. By using its *Soft Molding process* a rapid improvement in visual fields is obtained in the adaptation process.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	No
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	2.5 mm
Minimum VBOX	24 mm
Minimum fitting height (FH)	16 - 18 mm
Corridor	12 - 14 mm
Near reference point (NRP)	14 - 16 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 4.00 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

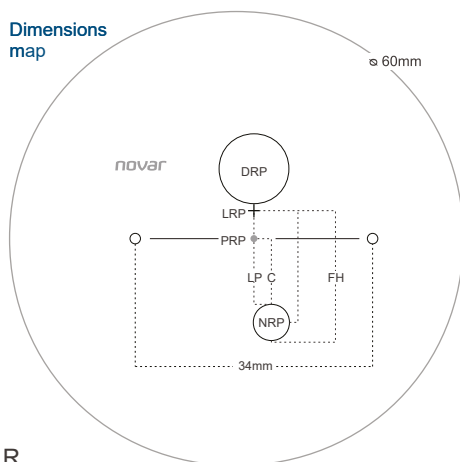


R

Power map

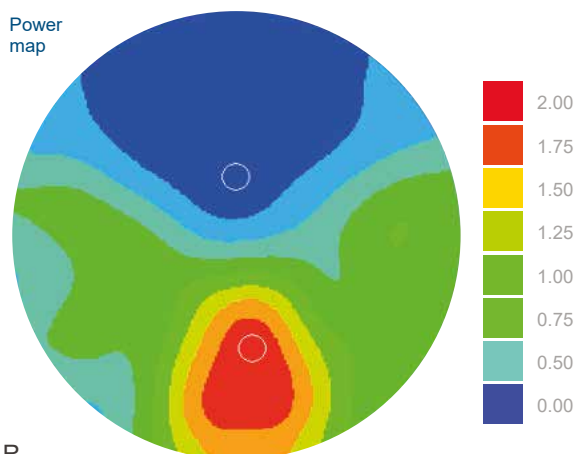

R

Cylinder map

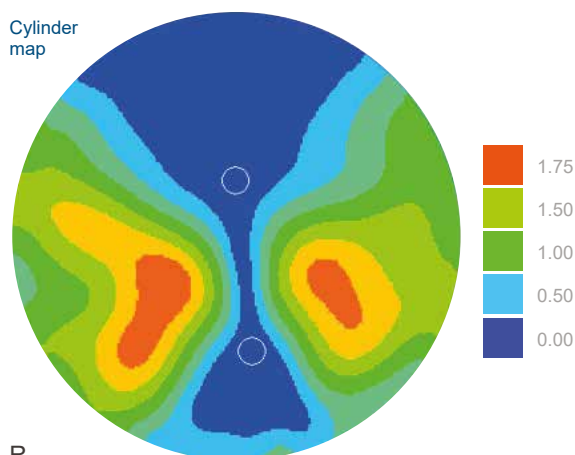

R

VISION



PRECISA PROGRESSIVE

General purpose progressive lens with balanced areas by using *Balanced Process* to achieve high performance at all distances.

Suitable for users who need lens for their daily tasks using all distances.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	Variable
Minimum VBOX	24 mm
Minimum fitting height (FH)	16 - 17 - 18 - 19 - 20 mm
Corridor	12 - 13 - 14 - 15 - 16 mm
Near reference point (NRP)	14 - 15 - 16 - 17 - 18 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 4.00 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes
Automatic corridor selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

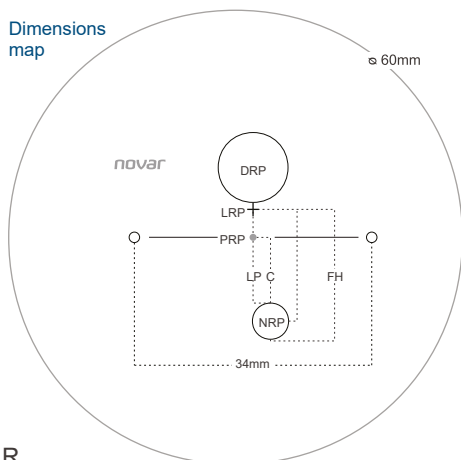


R

Power map

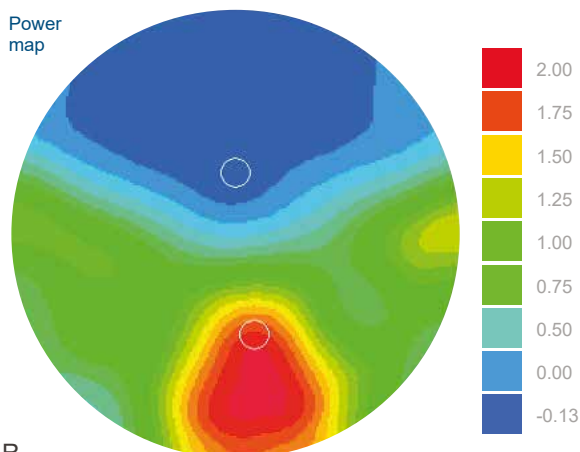

R

Cylinder map

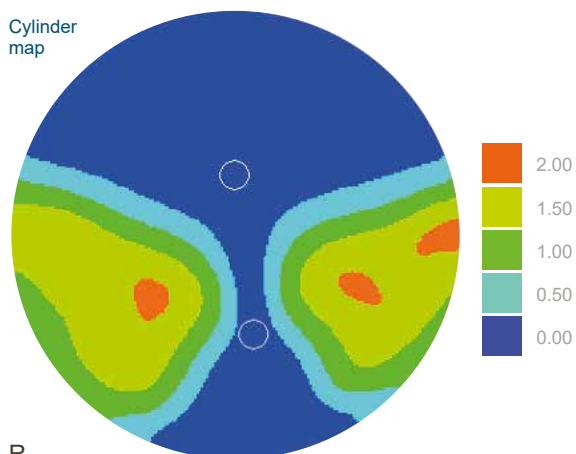

R

VISION



PRECISA SHORT PROGRESSIVE

Progressive lens with balanced design for small-sized frames.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Orgánico, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+6 mm
Layout reference point (LRP)	+2 mm
Inset	2 mm
Minimum VBOX	20 mm
Minimum fitting height (FH)	12 - 14 mm
Corridor	8 - 10 mm
Near reference point (NRP)	10 - 12 mm
Maximum diameter	75 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 4.00 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

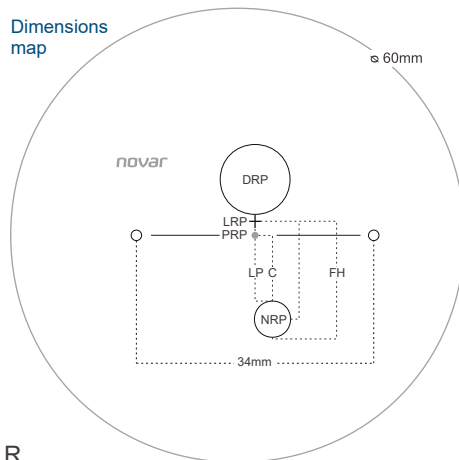


R

Power map

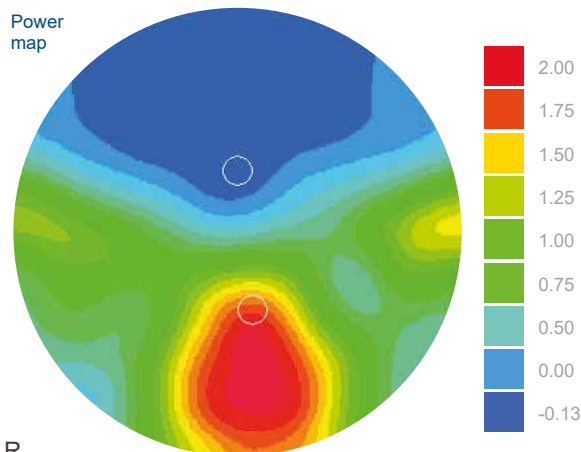

R

Cylinder map

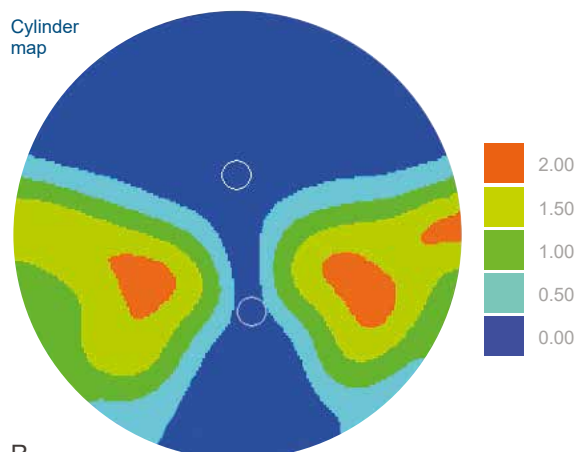

R

VISION





EVOLUTION PROGRESSIVE

Premium progressive lens designed with *Continuum Desing Technology & Smart Molding Process*. These technologies allow to place aberrations at the lower part of the lens to improve the breadth and to optimize all visual fields.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	Variable
Minimum VBOX	24 mm
Minimum fitting height (FH)	16 - 17 - 18 - 19 - 20 mm
Corridor	12 - 13 - 14 - 15 - 16 mm
Near reference point (NRP)	14 - 15 - 16 - 17 - 18 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.50 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes
Automatic corridor selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elipitical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

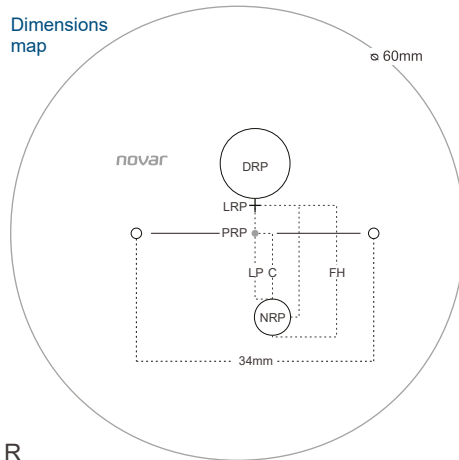


R

Power map

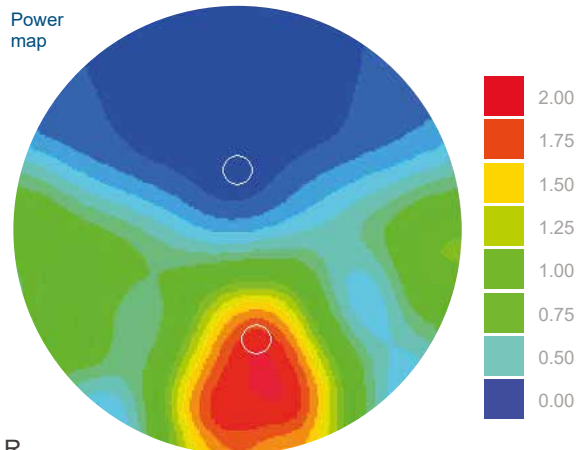

R

Cylinder map

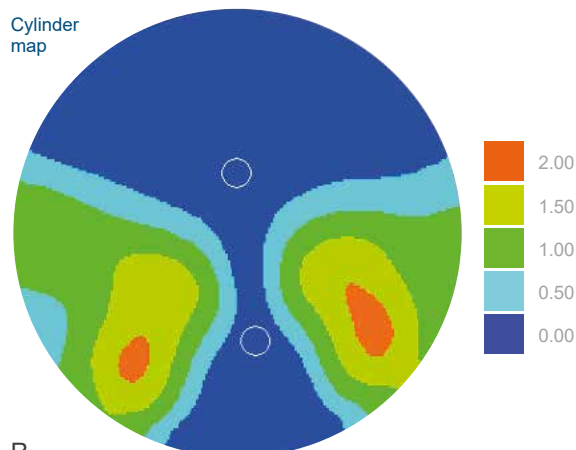

R


VISION



EVOLUTION SHORT PROGRESSIVE


The most evolved progressive lens for small frames.

Dimensions map

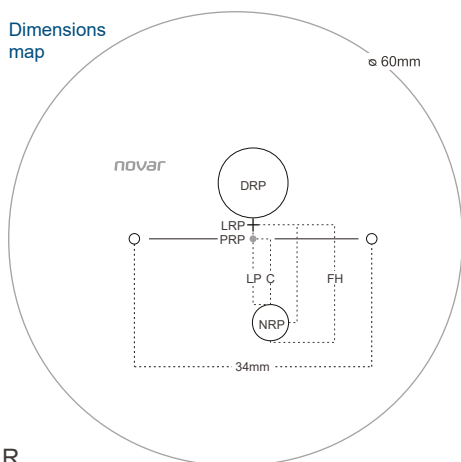


R

Power map

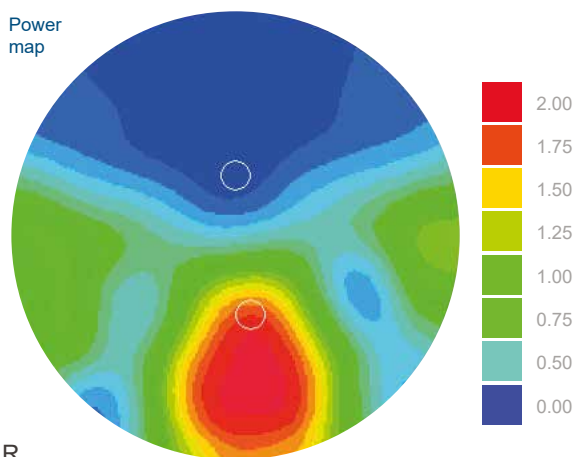

R

Cylinder map

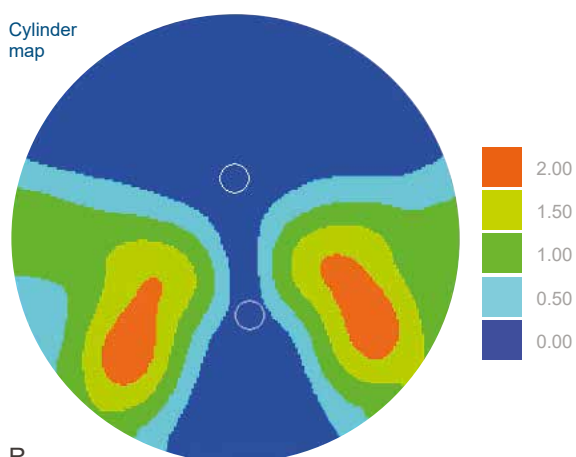

R

VISION



TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+6 mm
Layout reference point (LRP)	+2 mm
Inset	2 mm
Minimum VBOX	22 mm
Minimum fitting height (FH)	12-14 mm
Corridor	8 - 10 mm
Near reference point (NRP)	10 - 12 mm
Maximum diameter	75 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.50 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:

Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

eLIFE PROGRESSIVE

The *Adaptive Focus technology* improves the visual field widening the intermediate and near vision in order to face the challenges of modern life. Thus, a more comfortable and natural reading position is achieved by maintaining every feature of premium lenses in far vision.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	Variable
Minimum VBOX	22 mm
Minimum fitting height(FH)	16-17-18-19-20 mm
Corridor	12-13-14-15-16 mm
Near reference point (NRP)	14-15-16-17-18 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.50 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes
Automatic corridor selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

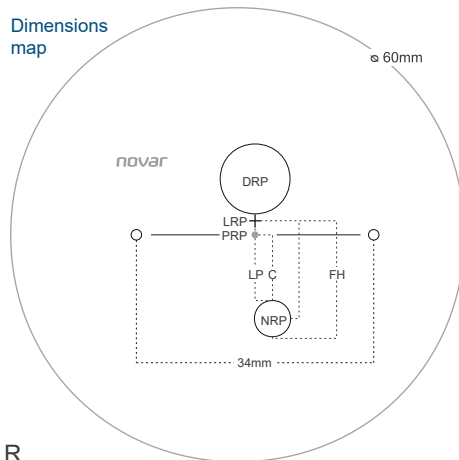


R

Power map

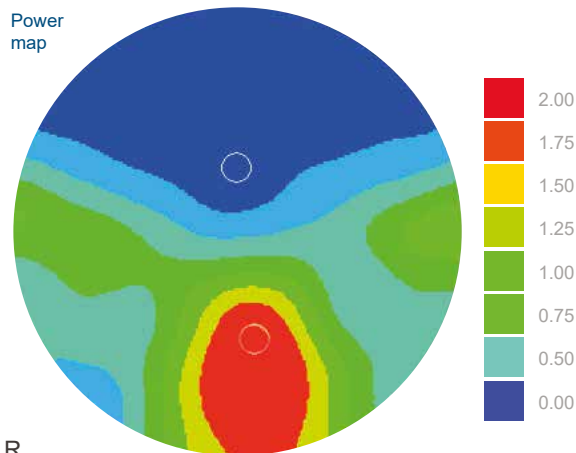

R

Cylinder map

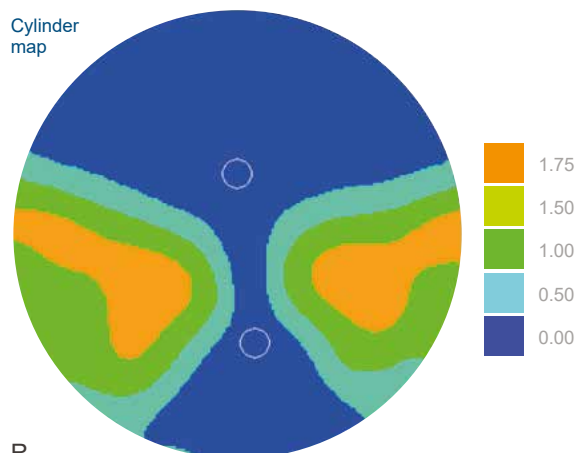


R

VISION




Dimensions map

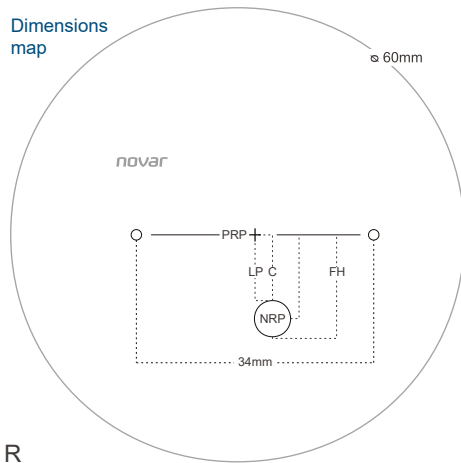


R

Power map

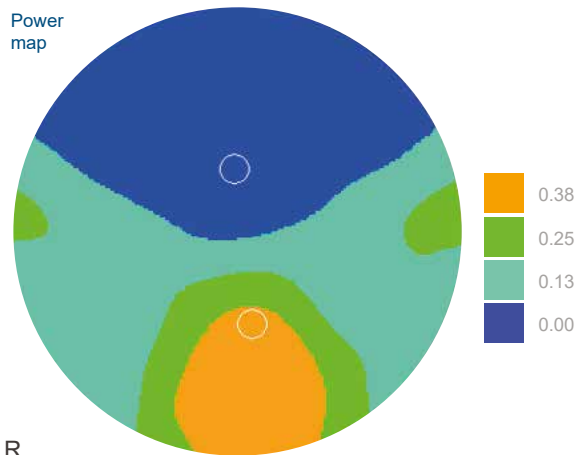


R

Cylinder map

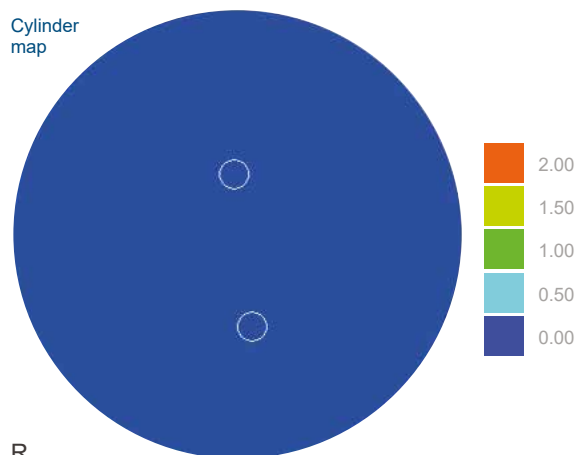

R

VISION

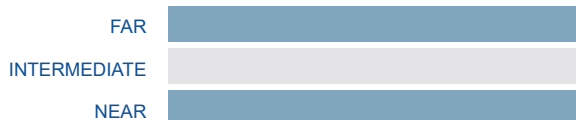

RELAX

Single vision lens developed for people aged 20-40 years who need to relax the eyes. The best lens for students and pre-presbyopic people who suffer from eye strain.

TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	0 mm
Layout reference point (LRP)	0 mm
Inset	2 mm
Minimum VBOX	16 mm
Minimum fitting height (FH)	16 mm
Near reference point (NRP)	10 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.36 / 0.52 / 0.72 / 0.96 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:

Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

DRIVE PROGRESSIVE

Progressive designed for those who spend most of their time behind the wheel. It incorporates *Free Periphery Process* which allows the lens to be free from peripheral astigmatism to achieve safer and more comfortable driving.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	Variable
Minimum VBOX	26 mm
Minimum fitting height (FH)	18 - 19 - 20 - 21 - 22 mm
Corridor	14 - 15 - 16 - 17 - 18 mm
Near reference point (NRP)	16 - 17 - 18 - 19 - 20 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.50 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

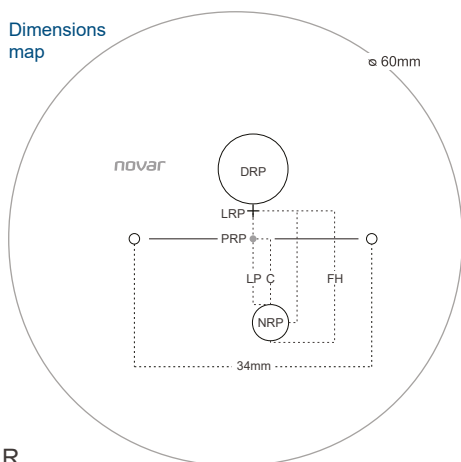


R

Power map

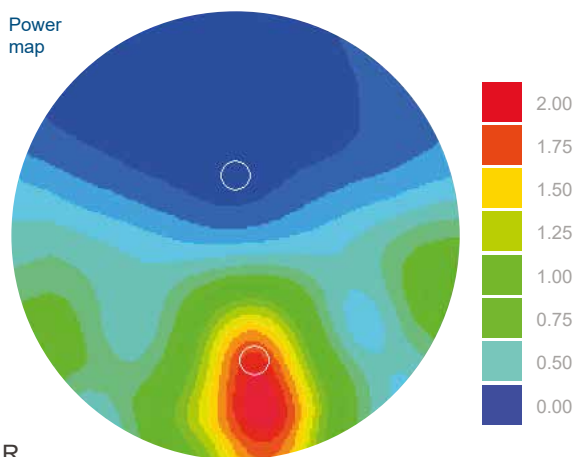

R

Cylinder map

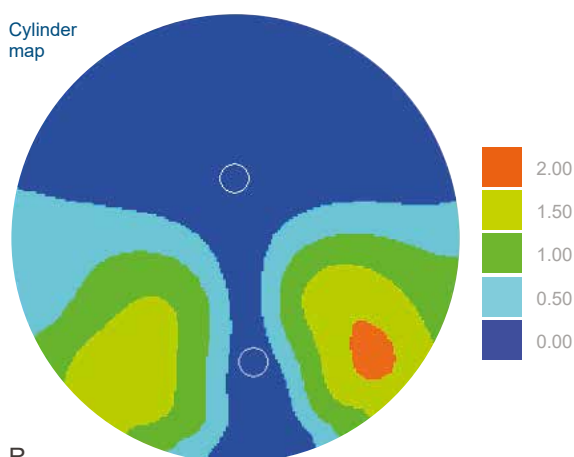

R

VISION



OUTDOOR PROGRESSIVE

Progressive designed for people who have a very active life outdoors.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	2.5 mm
Minimum VBOX	28 mm
Minimum fitting height (FH)	20 mm
Corridor	16 mm
Near reference point (NRP)	18 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.50 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

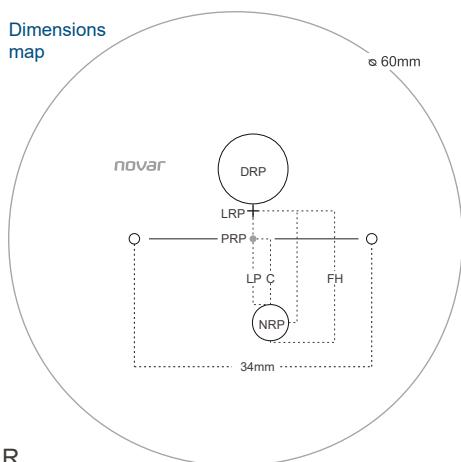


R

Power map

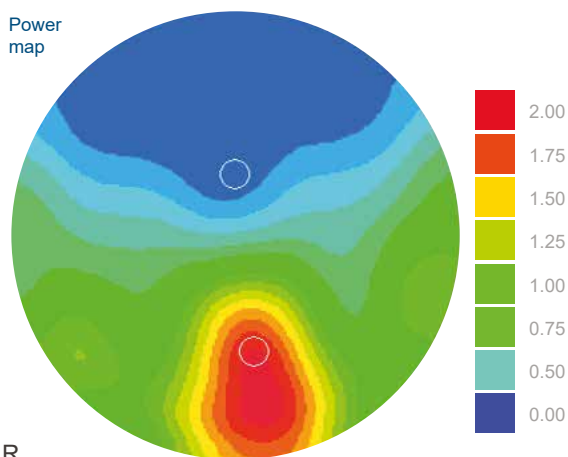

R

Cylinder map

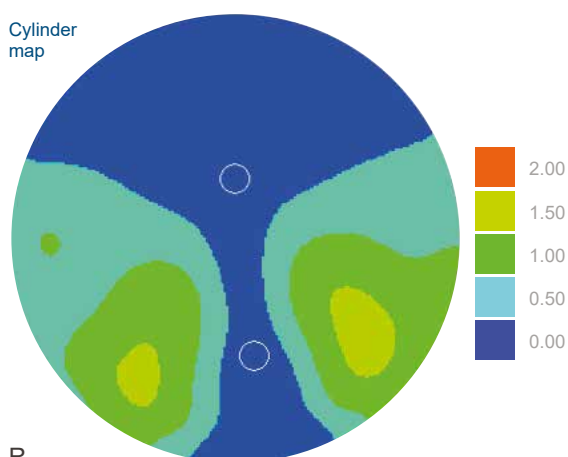

R

VISION



INDOOR PROGRESSIVE

Progressive lens designed for indoor work environments with great amplitude in near and intermediate vision.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	2.5 mm
Minimum VBOX	26 mm
Minimum fitting height (FH)	18 - 19 - 20 mm
Corridor	14 - 15 - 16 mm
Near reference point (NRP)	16 - 17 - 18 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.50 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

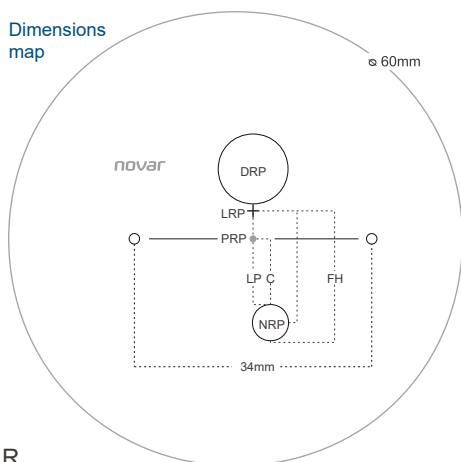


R

Power map

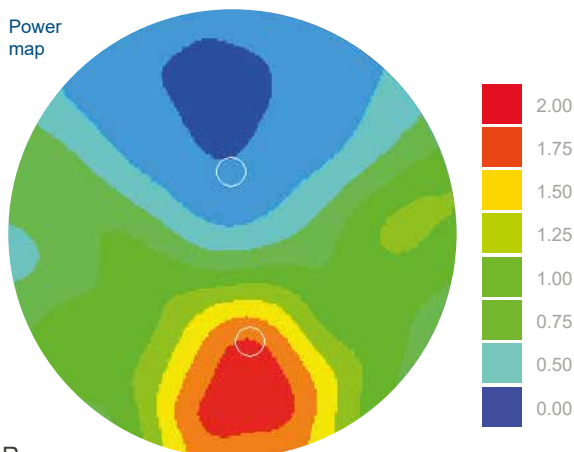

R

Cylinder map

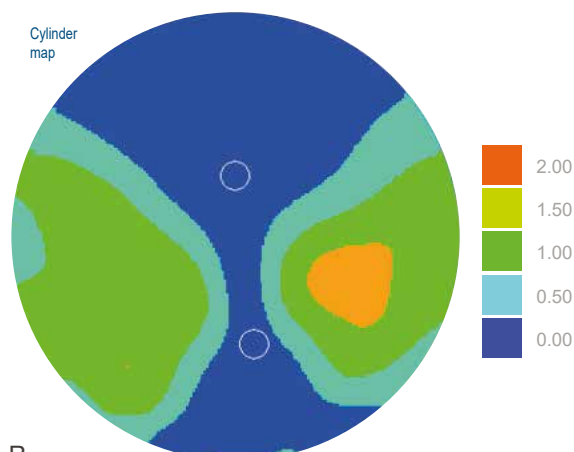

R

VISION



MONOVISION PROGRESSIVE

Progressive lens specially designed with null inset and specular symmetry for users with monocular vision due to convergence insufficiency or to the loss of an eye.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	Ø
Minimum VBOX	26 mm
Minimum fitting height (FH)	18 mm
Corridor	14 mm
Near reference point (NRP)	16 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.50 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

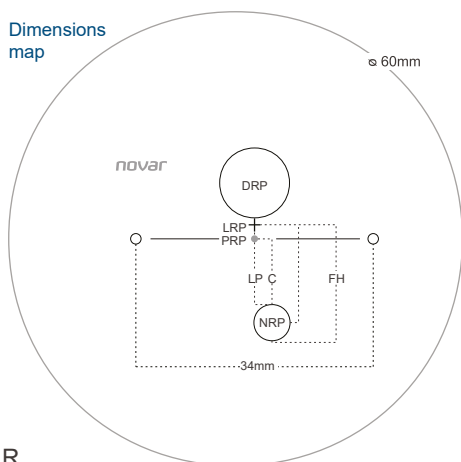


R

Power map

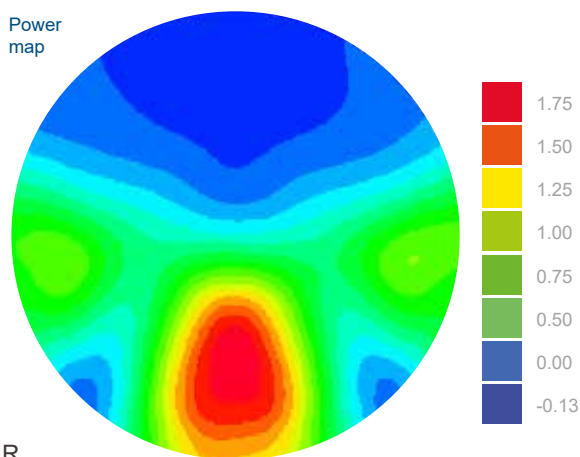

R

Cylinder map

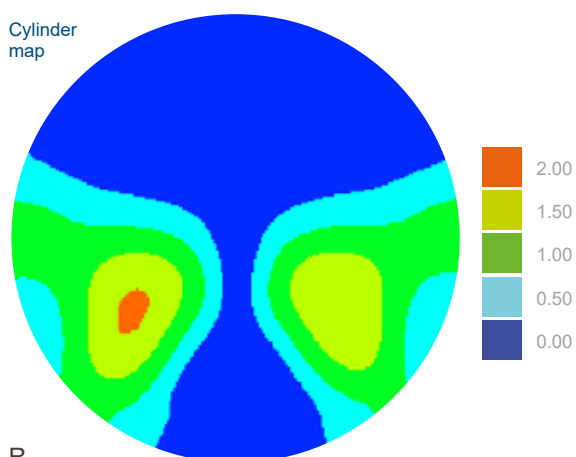

R

VISION



SENSITIVE PROGRESSIVE

Progressive lens designed with CDT & *Soft Molding* to achieve a soft design with an optimized periphery zone that meets the necessities of those who suffered from adaptation problems with other progressives due to their high distortions sensitivity.


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+8 mm
Layout reference point (LRP)	+4 mm
Inset	2,25
Minimum VBOX	22 mm
Minimum fitting height (FH)	18 mm
Corridor	14 mm
Near reference point (NRP)	16 mm
Maximum diameter	85 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	0.50 / 3.50 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes
Automatic corridor selection	Yes

Thickness calculation technology:


Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

Dimensions map

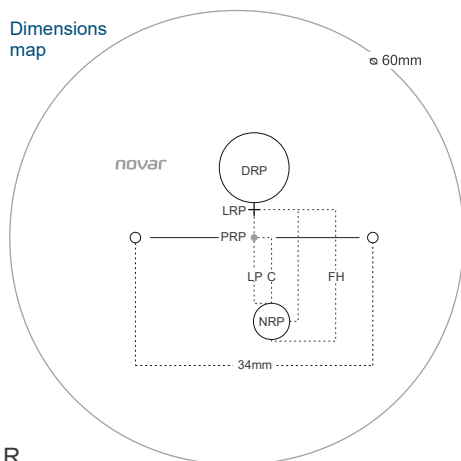


R

Power map

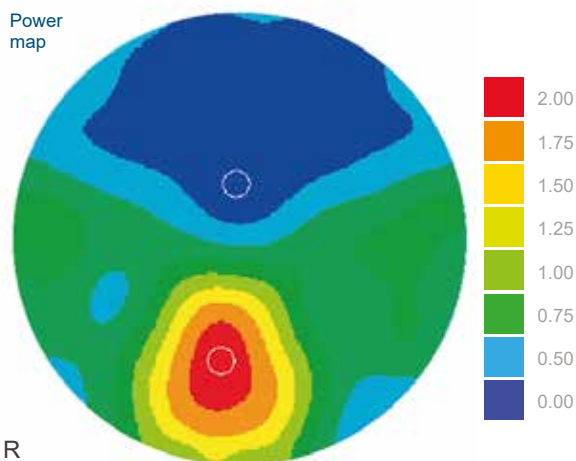

R

Cylinder map

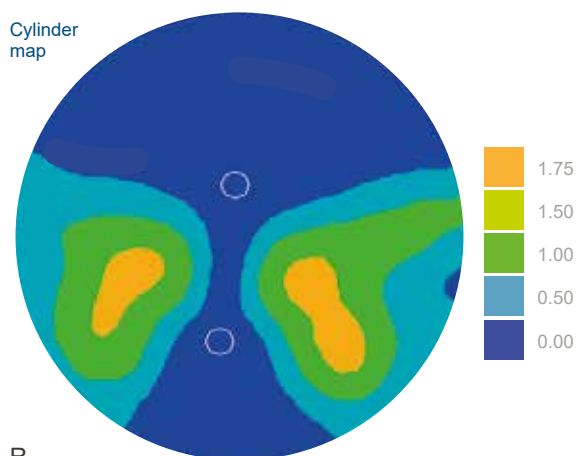


R

VISION


FAR	
INTERMEDIATE	
NEAR	

Dimensions map

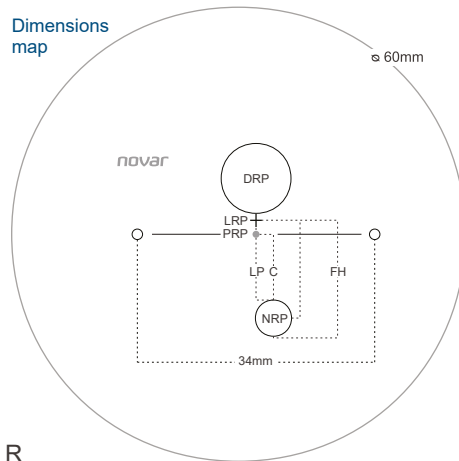


R

Power map

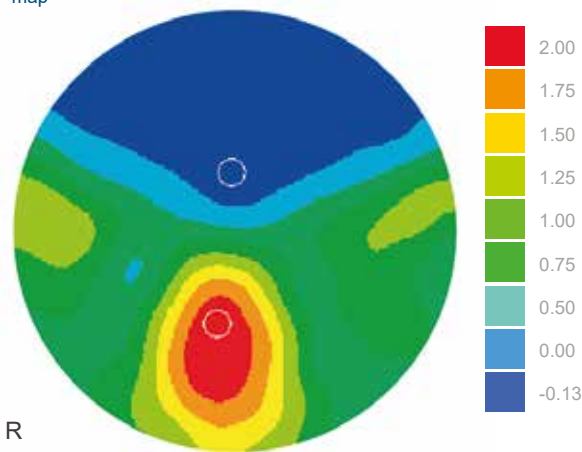

R

Cylinder map

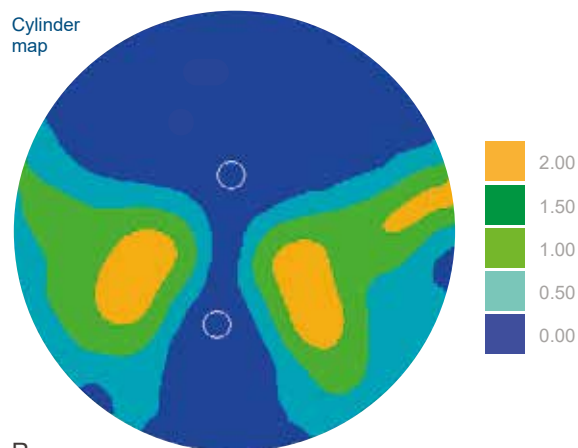

R

VISION



KIDCARE

Lenses designed to control the progression of nearsightedness in children


TECHNICAL DATA

Calculation technology	CDT® (Continuum Design Technology)
Surface resolution	4 K
Availability of materials	Organic, Trivex, Poli, High index
Availability of indexes	1.50, 1.53, 1.56, 1.59, 1.60, 1.67, 1.74
WEAR FIT customization	Yes
Distance reference point (DRP)	+6 mm
Layout reference point (LRP)	+2 mm
Inset	2 mm
Minimum VBOX	20 mm
Minimum fitting height (FH)	14 mm
Corridor	10 mm
Near reference point (NRP)	12 mm
Maximum diameter	75 mm
Spherical power range	-12 / +12 D
Cylindrical power range	-6 / +6 D
Addition range	2.00 D
Variable decentration	Yes
Prism prescribed by design	Yes
Automatic base curve selection	Yes

Thickness calculation technology:

Circular Fit	Yes
Elliptical Fit	Yes
Optimal Fit	Yes
Smart Fit	Yes
Blending Tech (Positives and Negatives)	Yes

SPECIFICATIONS & CHARACTERISTICS OF THE DESIGNS

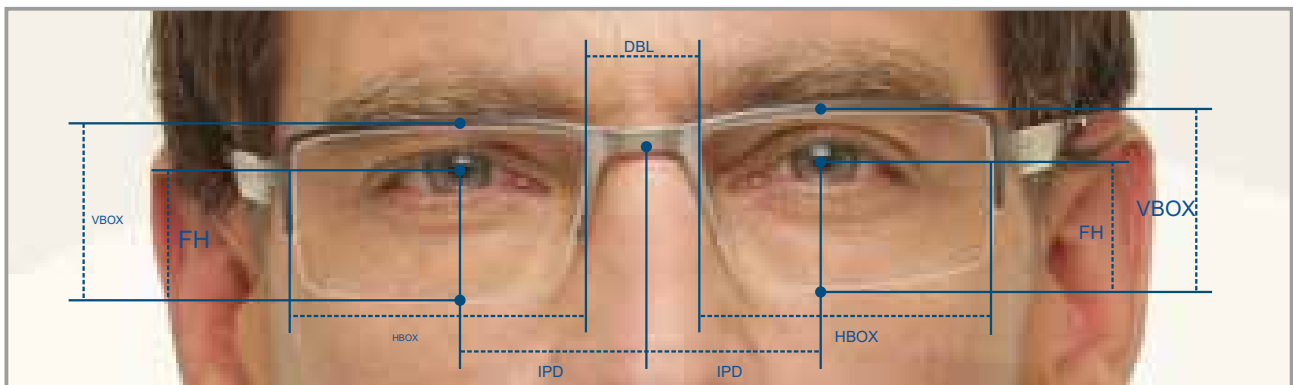

Calculation parameters

IPD - Interpupillary distance

SEGHT pupillary height

HBOX Horizontal BOX


DBL Distance Between Lenses - Bridge size

VBOX Vertical BOX

PANTO Pantoscopic tilt

ZTILT Wrapping angle

Technology for improving thickness

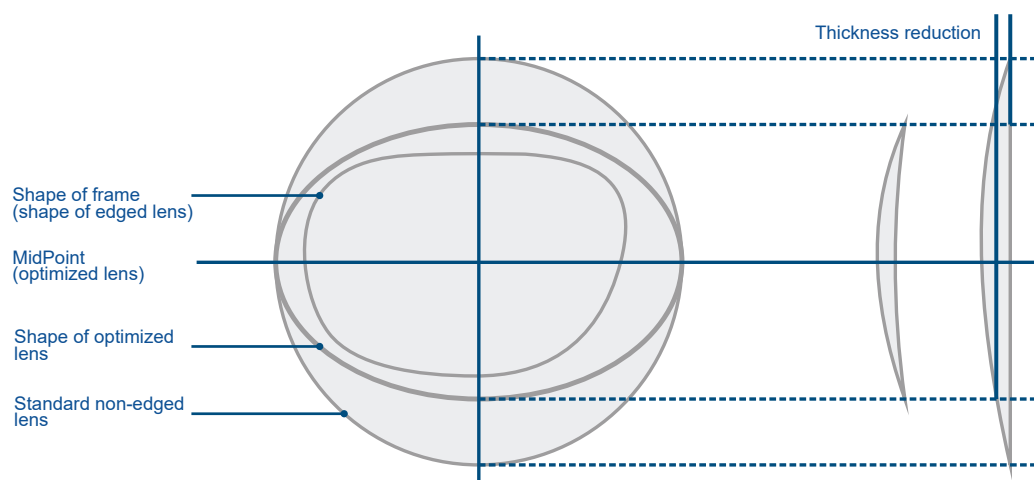

THICKNESS TECHNOLOGY

Calculation technology for the optimization and reduction of thickness for all kind of frames, shapes and materials.

CIRCULAR FIT

Thickness circular calculation from the lens diameter (CRIB) chosen by the user.

ELIPTICAL FIT

Optimization according to the frame horizontal width and vertical height so as to enhance cylindrical lenses thicknesses.

OPTIMAL FIT


Thickness reduction based on frame data (HBOX, VBOX, FED, FEDAX, DBL, IPD).

+ 25% Thickness reduction

SMART FIT

Smart thickness reduction based on the frame shape. (TRCFMT).

+ 33% Thickness reduction

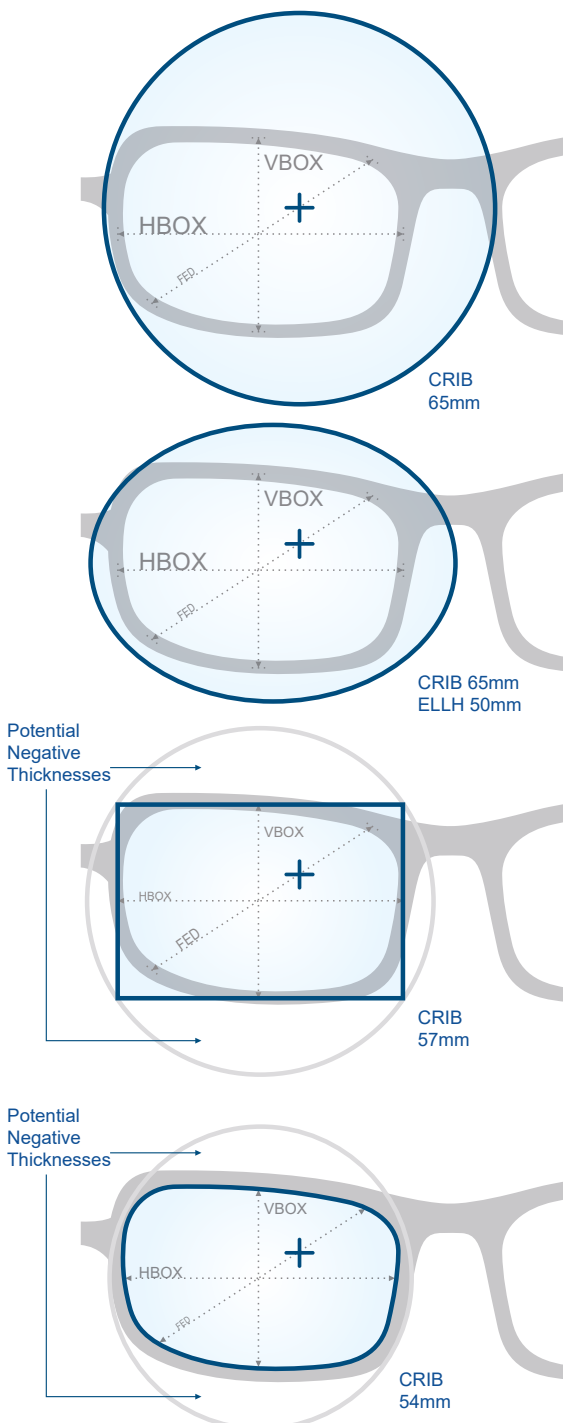


BLENDING TECHNOLOGY

LENTICULAR OPTICAL LENS

Technology developed to improve the edges of an ophtalmic lens. An optical zone is defined around the fitting cross and then, away from this zone, the thickness is significantly reduced with a gradual change of curvature in order to obtain a thinner lens.

This gradual change in curvature is available in STRONG or MODERATE version depending on the power of the lens and the frame shape.

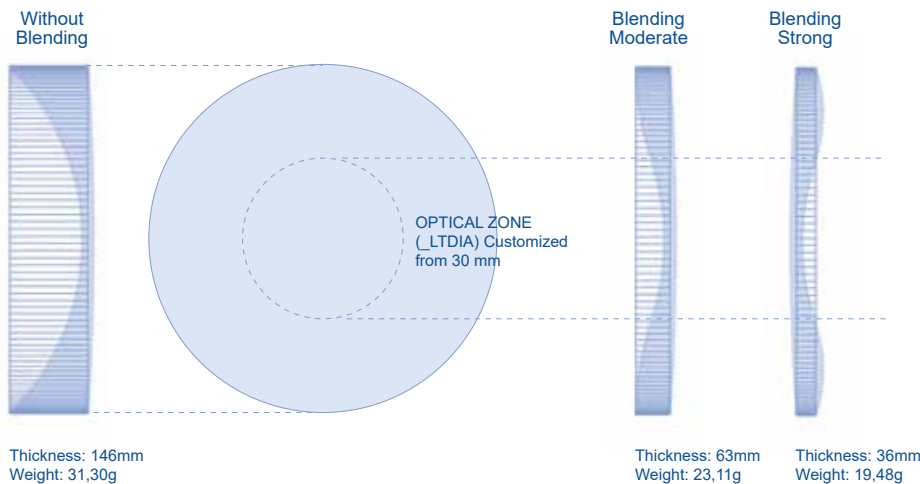

ADVANTAGES

- Thickness thinning
- More aesthetic lenses
- Available in versions: MODERATE and STRONG
- For positive and negatives lens
- For any design (even progressives lens)
- For any material
- Freedom to choose frames.

		EDGE THICKNESS			CENTER THICKNESS		
		material		%	material		%
		Poly 1.586	Org 1.499		Poly 1.586	Org 1.499	
		Esf.-6.00	Esf.-6.00		Esf.+6.00	Esf.+6.00	
		Base 0.50	Base 0.50		Base 9.00	Base 9.00	
NO BLENDING	_LTDIA	68mm	8.0	9.2	7.03	7.96	
STRONG	45mm	4.9	5.5	40	6.04	7.18	10
	40mm	4.3	4.8	48	5.56	6.69	17
	35mm	3.8	4.1	55	4.84	5.77	28
	30mm	3.3	3.5	62	4.06	4.81	39
MODERATE	45mm	5.1	5.7	38	6.35	7.53	5
	40mm	4.5	4.9	47	6.13	7.29	8
	35mm	3.9	4.3	53	5.84	6.97	12
	30mm	3.4	3.7	59	5.47	6.55	18


LENTICULAR FRAME SHAPE

Blending Frame (_LTFRAME) | Off (value 0), On (value 1)

Another lenticularization option is achieved by taking into account the frame shape. To achieve this, we define the millimetres to be applied in the zone surrounding the shape of the frame as explained in the graph.

Width and type: Available in MODERATE and STRONG version, the width of this zone is 5mm by default but it can be selected by the operator.

Blending Width (_LTWIDTH)
customized
from 1 to 10mm

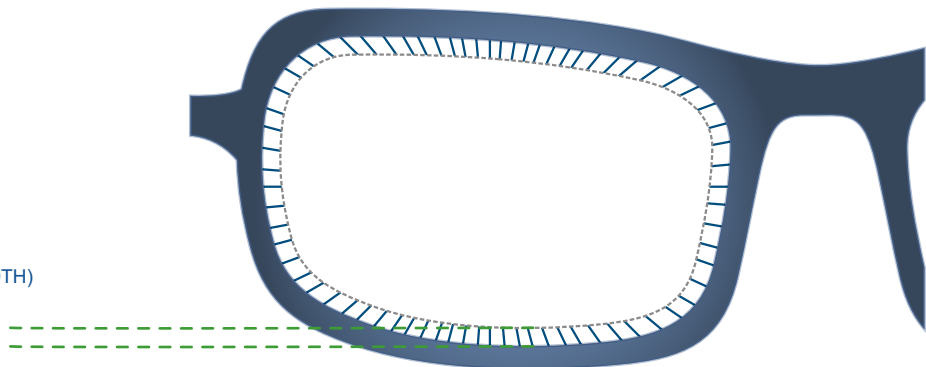



SOFTWARE DESIGNER



Software & hardware requirement specification (Hardware & Operating system)

Operating system (OS)

Windows XP, Windows Vista y Windows7/8, 32 y 64 bit, or Windows 2003

CPU

Intel or AMD (Processor speed 1.0 GHz or better)

RAM

minimum 1 GB

Processor

Inter Core 2 duo 2 GHz or similar

RAM

minimum 2 GB


HDD

100mb (minimum)

Internet access.

Standard VCA

NOVAR Designer follows the standard requirements of 3.08 for processing communications and file sharing "VCA Data Communication Standard v3.08".

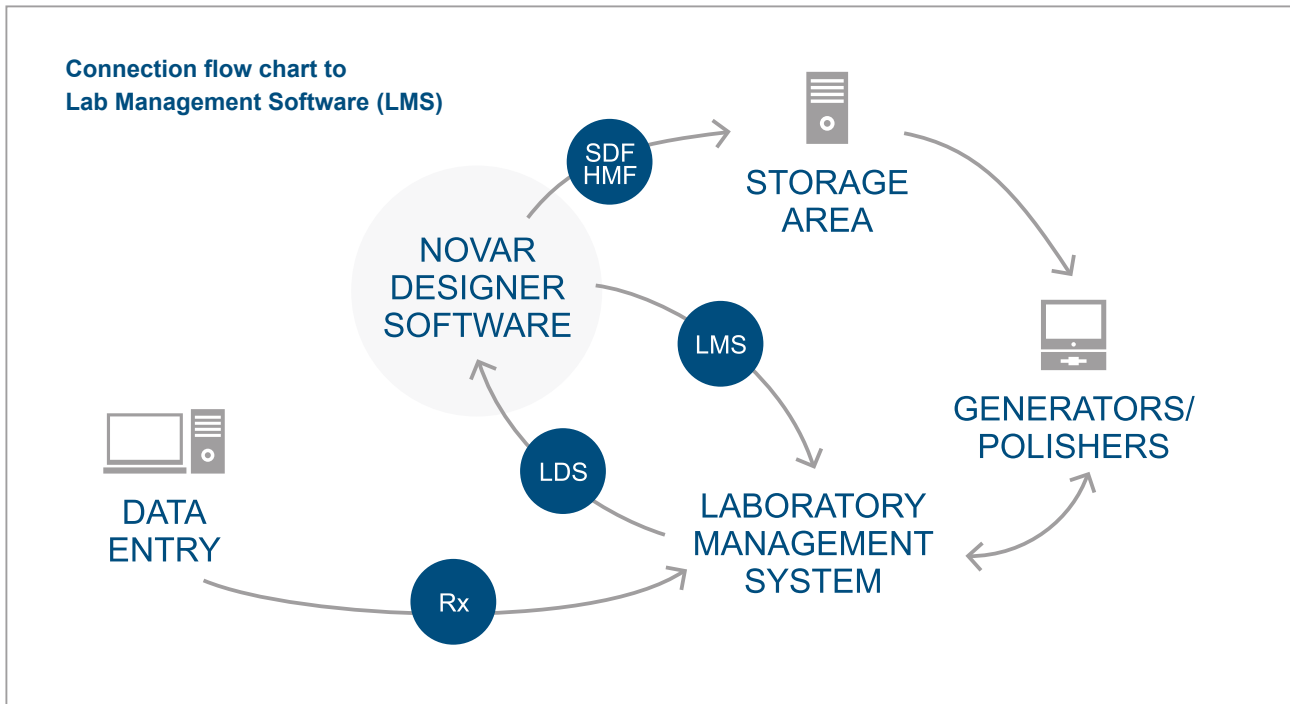

EQUIPEMENT & LMS

Designer developed under the AVC 3.08 standards.

Available for machines:

SATISLOH
COBURN
SCHNEIDER
OPTOTECH

With the following LMS:


RxUniverse
Innovations
Schneider LMS Basic / Plus
AfServer/Client
Rxdslab
SiuCalc
RxOffice y otros.

Download

Download the calculation software directly from our website and get the trial period:

www.novar-tech.com, FREE TEST section.

CALCULATOR VIEWER

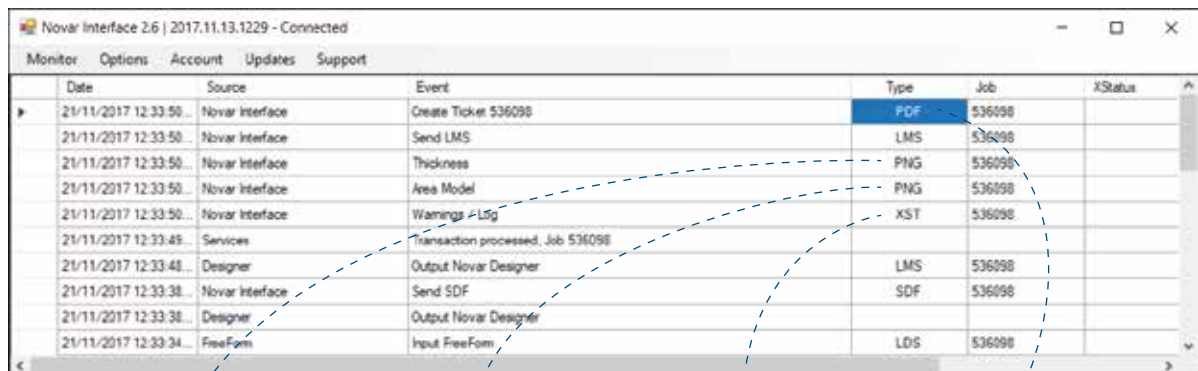

Date	Source	Event	Type	Job	Status
21/11/2017 12:33:50...	Novar Interface	Create Ticket 536098	PDF	536098	
21/11/2017 12:33:50...	Novar Interface	Send LMS	LMS	536098	
21/11/2017 12:33:50...	Novar Interface	Thickness	PNG	536098	
21/11/2017 12:33:50...	Novar Interface	Area Model	PNG	536098	
21/11/2017 12:33:50...	Novar Interface	Warnings - Lsg	XST	536098	
21/11/2017 12:33:45...	Services	Transaction processed, Job 536098			
21/11/2017 12:33:48...	Designer	Output Novar Designer	LMS	536098	
21/11/2017 12:33:38...	Novar Interface	Send SDF	SDF	536098	
21/11/2017 12:33:38...	Designer	Output Novar Designer			
21/11/2017 12:33:34...	FreeForm	Input FreeForm	LDS	536098	

A. THICKNESS
Graph

B. AREA MODEL
Graph

C. Warning
& Notifications


D. Calculation Overall
Report

REPORTS FOR EACH CALCULATION

A. THICKNESS graph


Front view of the lens with all its thicknesses (nasal, temporal, superior, inferior, minimum -in green- and maximum -in red-).

3D profile view with temporal thickness and cross section. In addition, final diameter information is displayed together with decentration and weight data (in grams) of the lens. Graphs available in test calculation. The goal is to help the lab operator to get more information when it comes to improving the thicknesses of the lens to be designed.

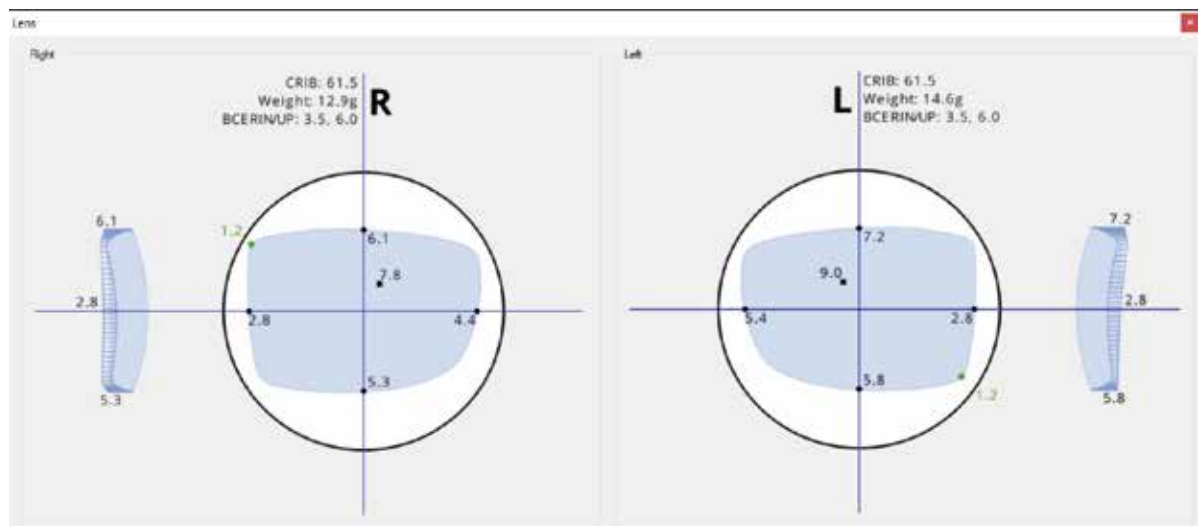

B. AREA MODEL graph

It allows to establish the exact place where the far and near vision areas will be located within the frame.

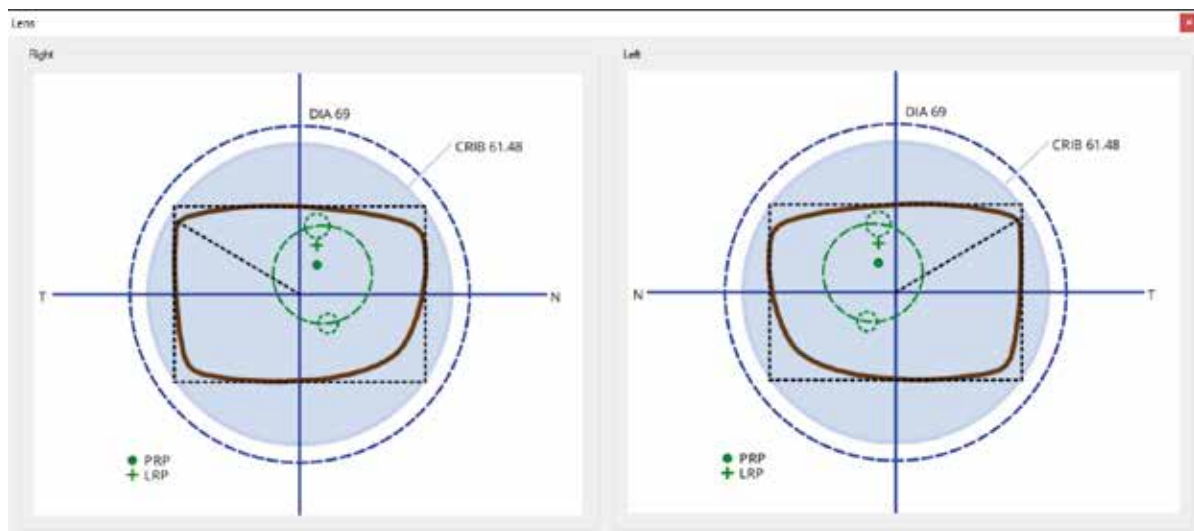

C. Warning & Notifications


Each calculation displays information about different topics related to the calculation such as: ideal bases, prisms, blending, etc.



D. Calculation overall report

Report to be delivered to the opticians-optometrists along with their lenses with information related to the original prescription, calculated prescription and its adjustments, in case of using the Wear Fit customization option.

In the center of the report, the customized centration chart is printed with information on thickness, decentration and weight on its sides. At the bottom, customized maps of power, astigmatism and thickness of the lens will be displayed.

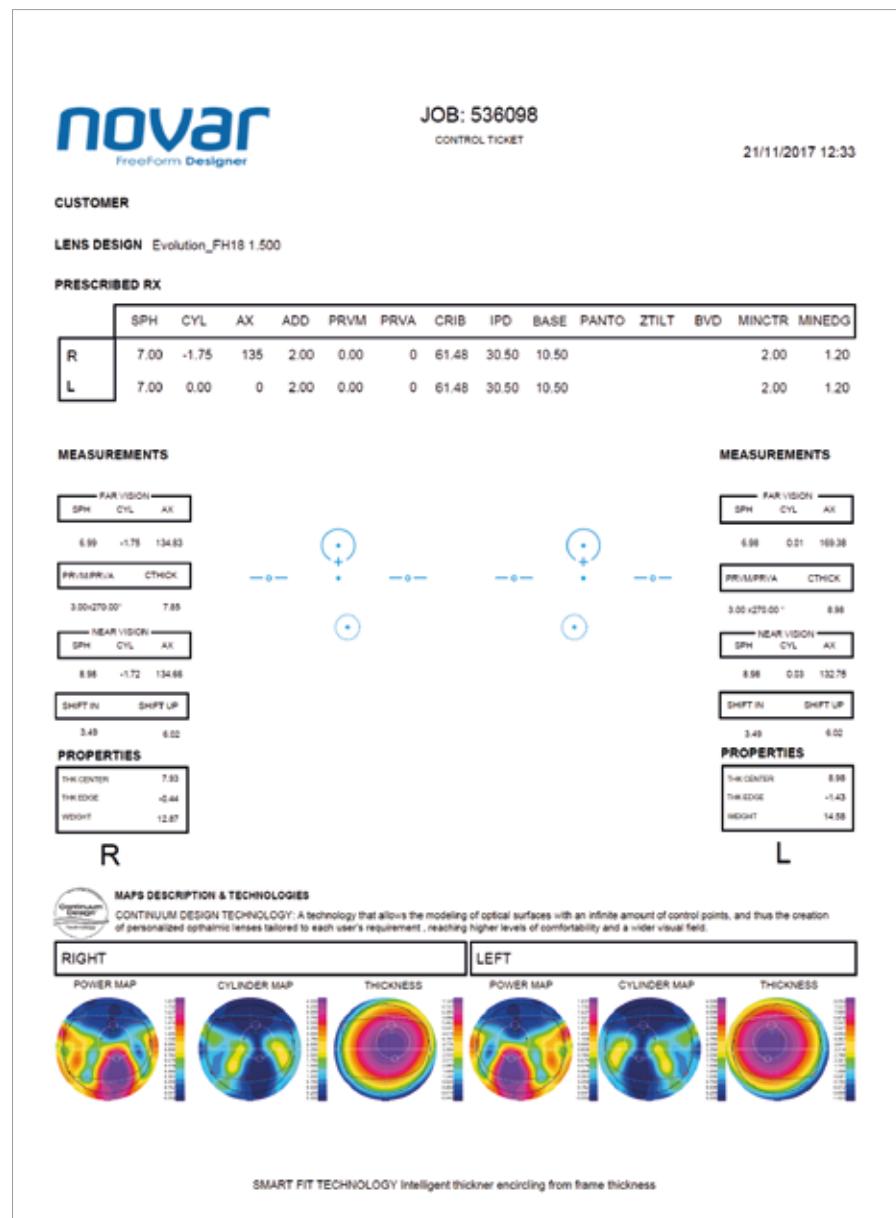


MARKETING


Graphic material for advertisements: posters, banners, flyers, marquees, etc - Advertising videos - points of sale material. All Communication Materials to be downloaded.

www.novar-tech.com, Marketing section.

Mkt
TV/Radio spot
Commercials

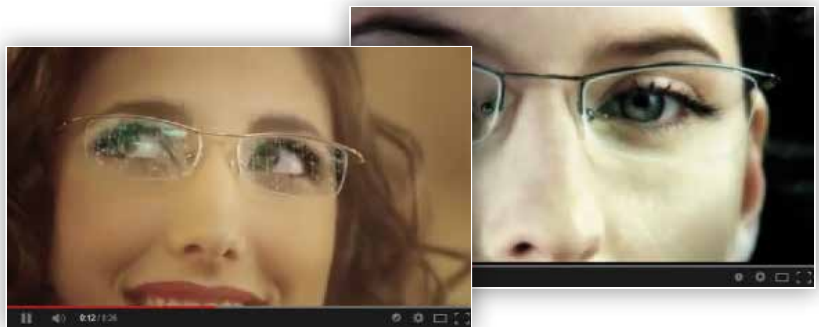

Mkt
Videos

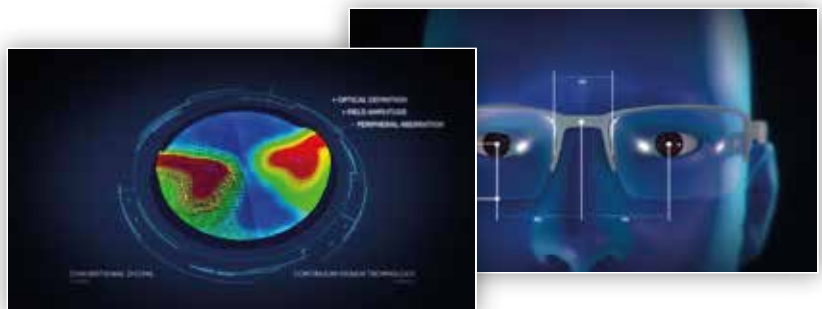

Mkt
Display
Advertising
For Opticians



Mkt Display Advertising For consumers



Mkt POP



+ INFO
www.novar-tech.com



DIGITAL CENTER

Accurate measure of IPD.

With *Mirroring Alignment System*.

Real thickness calculation.

Real time connection to the Freeform lens

Designer software to get information in advance.

MKT Tools.

Demonstrations about designs, technologies and materials for clear explanations.

